


# Optimal Solett

Preventative mattress replacement for  
lightweight and sensitive patients

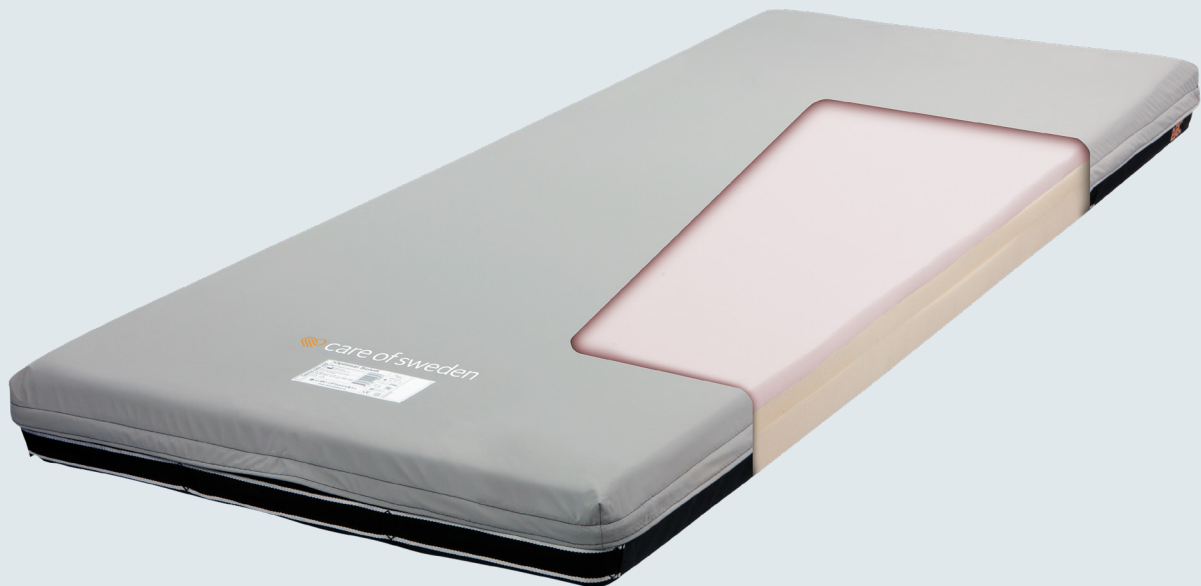

Optimal Solett with Olivia/Evac cover

# Optimal Solett

Optimal Solett is used as an aid in the prevention and treatment of pressure ulcers up to category 1<sup>(1)</sup>. It has been specially designed for lightweight and sensitive users. This mattress replacement provides good pressure redistribution.

Optimal Solett is delivered with the Olivia/Evac cover. It is equipped with a divisible, drip shield zip and feature handles for evacuation.

Care of Sweden's mattresses are supplied with a four way stretch hygienic outer cover which conforms to relevant standards for hygiene, quality and safety<sup>(2)(4)</sup>. All covers are manufactured from modern, recyclable materials and are vapour permeable<sup>(3)</sup>, which reduce the risk of skin maceration.

## SIZES

80/85/90/105/120x200 cm and 80/85/90x210 cm.  
If you need any other sizes please contact your distributor.



## AVAILABLE COVERS

### Olivia/Evac Cover

Divisible, drip shield zip and handles for evacuation on one side.

### Cover Olivia

Drip shield zip on three sides.

### Cover Nordic

Welded seams, drip shield zip on three sides and handles on one side.

### Steel/Evac Cover

2-piece cover with handles for evacuation. Welded seams and a protected zipper on all sides.

For need of non hygienic cover, please contact your distributor.


## REFERENCES

(1) European Pressure Ulcer Advisory Panel, National Pressure Injury Advisory Panel and Pan Pacific Pressure Injury Alliance. Prevention and Treatment of Pressure Ulcers/Injuries: Clinical Practice Guideline. The International Guideline. Emily Haesler (Ed). EPUAP/NPIAP/PPPIA: 2019

(2) CDC (2003), Guidelines for Environmental Infection Control in Health-Care facilities, updated 2019

(3) SS-EN ISO 15496:2004, DIN 53122-1

(4) SS 876 00 20:2017, Healthcare textiles - Mattresses - Specifications and requirements

FEATURES & CHARACTERISTICS	FACTS
<p><b>Pressure</b> Elastic material with pressure equalisation that provides good comfort.</p> <p><b>Shear forces</b> The mattress is specially designed to minimise the effects of shear forces.</p> <p><b>Moisture</b> The hygienic cover is vapour permeable<sup>(3)</sup>, which reduces the risk of skin maceration.</p> <p><b>Hygiene</b> The cover is easy to clean and has a good documented hygienic standard<sup>(2)(4)</sup>.</p>	<ul style="list-style-type: none"><li>• Core: HR-foam 35 kg/m<sup>3</sup> + visco-elastic foam 50 kg/m<sup>3</sup></li><li>• Height: 12 cm</li><li>• Function: prevention and treatment up to and including category 1<sup>(1)</sup></li><li>• Recommended user weight: up to 60 kg</li><li>• Fire: EN 597-1, -2</li><li>• Warranty: 2 years</li><li>• Cleaning of cover: wipe with cleaning agent and/or disinfectants</li><li>• Machine wash max 95 °C, tumble drying</li></ul> <p> Always read the instruction for use prior to usage of the product.</p> <p>Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2023. All rights reserved. Last update: 2023-09-22.</p>