


# CuroCell® AUTO4 CX13

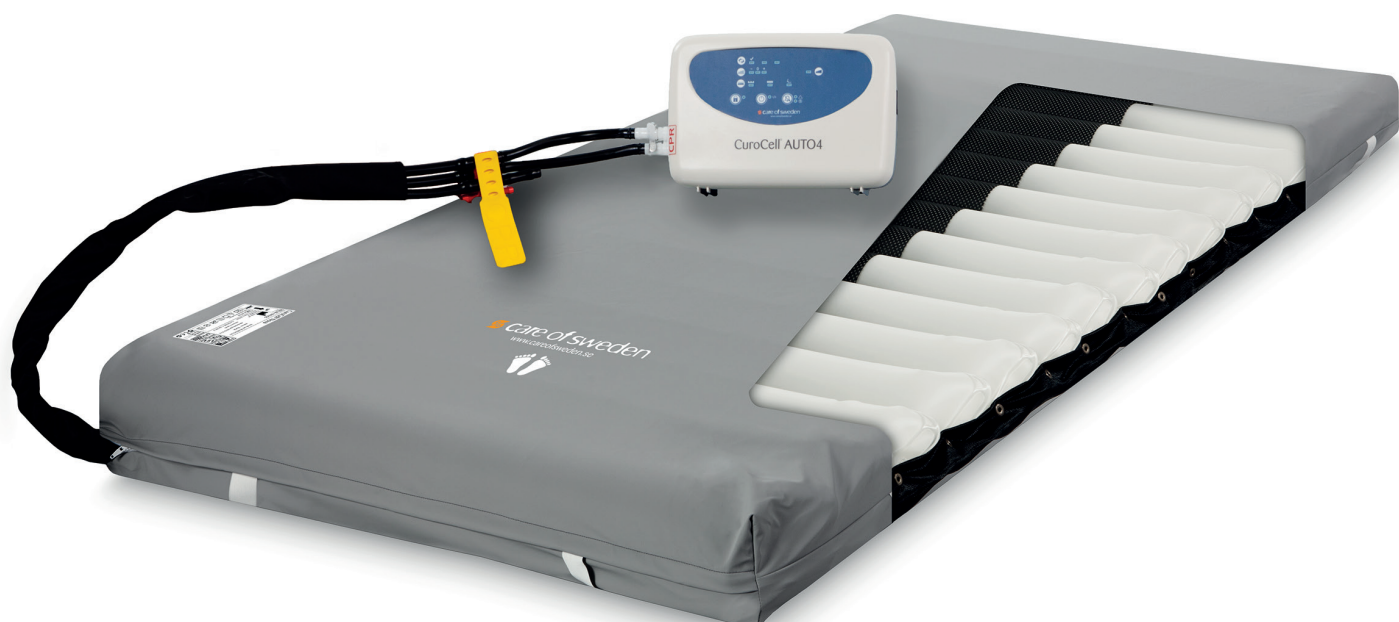

Effective and functional mattress overlay

# CuroCell® AUTO4 CX13

CuroCell® AUTO4 CX13 is used preventively and as an aid in the treatment of pressure ulcers up to category 4<sup>(1)</sup>. It can be used in all types of medical environments and is a highly-functional and cost-effective mattress system.

CuroCell® AUTO4 CX13 can be used with alternating pressure or static low pressure. The pump, which is easily maintained and managed, is also equipped with an automatic function that adjusts the mattress to the user's weight, posture and position.

The mattress has replaceable cells that are fitted inside a separate inner cover.

Care of Sweden's mattresses are supplied with a four way stretch hygienic outer cover which conforms to relevant standards for hygiene<sup>(2)</sup>, quality and safety. All covers are manufactured from modern, recyclable materials and are vapour permeable<sup>(3)</sup>, which reduce the risk of skin maceration.

## Facts

Type: Alternating pressure mattress

Height: 13 cm

Function: Prevention and treatment up to category 4<sup>(1)</sup>

Recommended user weight: 0-180 kg

Fire: EN 597-1, -2

Power supply: 230 V/50 Hz

Sound level, pump: 21 dBA

Warranty: 2 years

Cleaning covers: Wipe clean with detergent and/or disinfectant.

Machine wash max. 95°C, tumble dry.



## Functions & Characteristics

- CPR function – rapid deflation in case of CPR
- Digital pump with self-adjusting technology
- Easy to use
- Handles for transport and moving
- Audio/visual low pressure alarm
- Audio/visual power loss alarm
- Soft, elastic hygiene cover that is vapour permeable<sup>(2)</sup>
- Static function with automatic return
- Transport function
- Quiet and energy-efficient pump
- Replaceable polyurethane air cells

## Sizes

80/85/90/105/120x200 cm and 80/85/90x210 cm.

If other sizes are required, please contact Care of Sweden.

## Options

- Heel function
- Mesh fabric – for additional comfort
- Hygiene cover with welded seams
- Cable holder
- Sturdy transport bag
- Sheet holder
- Low friction cover

## References

(1) National Pressure Ulcer Advisory Panel, European Pressure Ulcer Advisory Panel and Pan Pacific Pressure Injury Alliance. Prevention and Treatment of Pressure Ulcers: Quick Reference Guide. Emily Haesler (Ed.). Cambridge Media: Osborne Park, Western Australia; 2014.

(2) ISO 16603 "resistance to penetration by blood and body fluids", ISO 16604 "resistance to penetration by blood-borne pathogens", EN 14126, SS 876 00 19 "Health care textiles – Bacteria penetration – Wet".

(3) SS-EN ISO 15496:2004, DIN 53.122.



Read the user manual prior to using the product.

Care of Sweden and CuroCell are registered trademarks of Care of Sweden AB.  
© Care of Sweden AB, 2016. All rights reserved. Last updated: 2016-11-23.



### Mattress selector

– Guidelines for selecting a product

Read about how to choose the right mattress at:

[www.careofsweden.com/mattress-selector](http://www.careofsweden.com/mattress-selector)

 care of sweden  
Supporting Life

Care of Sweden AB, Box 146, SE-514 23 Tranemo, Sweden  
+46 (0)771 106 600 | [info@careofsweden.se](mailto:info@careofsweden.se) | [www.careofsweden.com](http://www.careofsweden.com)