

Curera[®] Seat Cushion

Stable & comfortable seat cushion




Curera® Seat Cushion

Curera® Seat Cushion is a stable and comfortable cushion for wheel chairs, but can also be used on furniture. It has a cube design, and the many, small cubes are independently agile and conform to the shape of the body.

The cubes provide support and equalize the high pressure over the entire sitting surface. Tenderness and pain is alleviated and pressure ulcers are prevented.

The washable core is produced from foam with high flexibility and density, in order to make it stable and durable. The cushion has a welded and flexible hygiene cover, which has a drip-shielded zip fastener on the rear side. The hygienic cover is detachable,

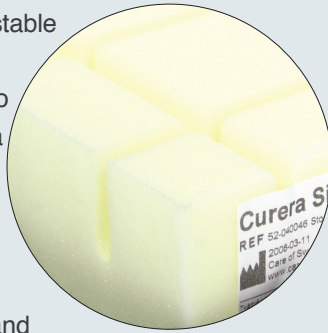

liquid resistant, vapour permeable and breathes, which prolongs the lifespan of the cushion. The cover is easily wiped off or machine washed and it reduces the need for washing the cushion itself.

Thanks to the cube design, the cushion can be used in a hanging seat. However; if the seat hangs substantially, it is recommended to use an equalization board.



Always read the instruction for use prior to usage of the product.

Care of Sweden are trademarks of Care of Sweden AB. © Care of Sweden AB, 2021. All rights reserved. Last update: 2021-06-04.

FEATURES & CHARACTERISTICS	FACTS
<ul style="list-style-type: none"> • Cube design • Easy cleaning • Unique labelling • Welded cover 	<ul style="list-style-type: none"> • Filling: 50kg/m3 polyether foam • Cover: Polyurethane covered polyester • Height: 6 cm • Fire: ISO 7176-16 • Warranty: 1 year • Classification: ISO code 18 10 06 (pressure redistributing) • Cleaning of cover: Wipe with cleaning agent and/or disinfectants. • Machine wash max 95 °C, tumble drying.
<p>SIZE</p>	
<p>35x35x6, 35x41x6, 40x41x6, 40x46x6, 40x50x6, 45x41x6, 45x46x6, 45x50x6, 50x46x6, 50x50x6 cm.</p>	