

Curera[®] Contur Seat Cushion

Comfortable and popular seat cushion



Curera[®] Contur Seat Cushion with Standard Cover

Curera® Contur Seat Cushion

Curera® Contur is a comfortable and pressure redistribution seat cushion made from an anatomically moulded base of polyurethane foam with memory foam in the seat well

The anatomical design distributes the pressure away from the underside of the thighs, provides a good and stable sitting position, good pressure redistribution for the seat bones and reduced risk of shearing. Recesses for lumbar supports, rounded rear edge and leg divider also give it good stability and fit in a wheelchair. Contur provides a good level of sitting comfort and is suitable for people sitting for both longer and shorter periods.

Contur has a waterproof surface that prevents liquids from passing through and is easy to clean. The cushion comes with a durable cover in grey polyester plush with high friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.



Always read the instructions for use prior to use.

Care of Sweden is a trademark of Care of Sweden AB.
© Care of Sweden AB, 2021. All rights reserved. Last update: 2021-06-03.

AVAILABLE COVERS

Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

Deluxe Cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

Onyx Cover

Material made of 4-way stretch on the top side and an anti-slip underside. Comes with handles on the underside. Available with a mesh layer for a better microclimate and a generally higher comfort.

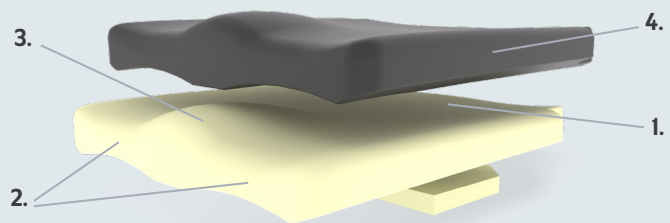

1. SEAT WELL

2. RECESSES FOR LEGS

3. LEG DIVIDER

4. WATERPROOF SURFACE

PROPERTIES	FACTS
<ul style="list-style-type: none"> Anatomically moulded Easy to clean Accessories: Levelling plate 	<ul style="list-style-type: none"> Height: 5 cm Warranty: 2 years Material in the cushion: Polyurethane foam, memory foam, vulcanised rubber Material in the cover: Polyester/polyurethane Classification: ISO-code 18 10 06 - Pressure redistribution Fire: Complies with EN 1021-1, EN 1021-2, ISO 7176-16

PASSFORMSGUIDE:

(Recommended max user weights in KG and available standard sizes. Other sizes than listed below can be ordered. Please contact your distributor.)

Width cm		30	33	35	38	40	42	45	47	50	52
Depth cm	30	70	70								
	35	70	70	70		100					
	40	70	70	70	100	125	125	125	125	150	150
	42					125	125	125	125	150	
	45			70	100	125	125	125	150	150	150
	47					125	125	125			
	50				100	125	125	125		150	150
	55									150	