

Curera[®] Plantur Seat Cushion

Seat Cushion with flat top side



Curera[®] Plantur Seat Cushion with Standard Cover

Curera® Plantur Seat Cushion

Curera® Plantur is a pressure redistribution seat, made from moulded polyurethane foam with a flat top side and with memory foam filling on the underside.

The seat well and recesses for legs distribute the pressure away from the underside of the thighs, providing a good sitting position, good pressure redistribution for the seat bones and reduced risk of shearing. Recesses for lumbar supports and rounded rear edge also give it good stability and fit in a wheelchair. Plantur provides comfort suitable for people sitting for both longer and shorter periods.

Curera® Plantur has a waterproof surface that prevents liquids from passing through and is easy to clean. The cushion comes with a durable cover in grey polyester plush with high friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.



Always read the instructions for use prior to use.

Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2021. All rights reserved. Last update: 2021-06-07.

AVAILABLE COVERS

Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

Onyx Cover

Material made of 4-way stretch on the top side and an anti-slip underside. Comes with handles on the underside. Available with a mesh layer for a better microclimate and a generally higher comfort.

Onyx Cover+

Flexible, elastic and removable cover made in 4-way stretch on the top and anti-slip on the lower part which counteracts forward sliding. Intermediate layer of mesh which contributes to a good microclimate and comfort. Handle on the underside.

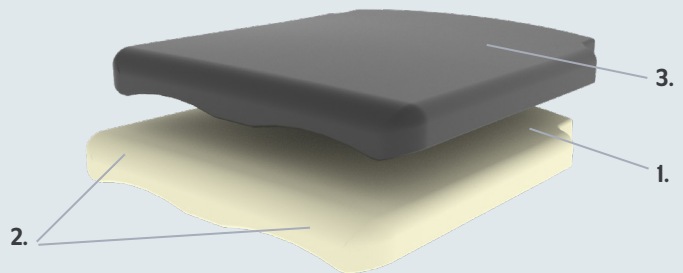

1. SEAT WELL
2. RECESSES FOR LEGS
3. WATERPROOF SURFACE

PROPERTIES	FACTS
<ul style="list-style-type: none"> • Easy to clean • Flat surface 	<ul style="list-style-type: none"> • Height: 5 cm • Warranty: 2 years • Material in the cushion: Polyurethane foam, memory foam, vulcanised rubber • Material in the cover: Polyester/polyurethane • Classification: ISO-kod 18 10 06 - Pressure redistribution • Fire: SS-EN 1021-1, SS-EN 1021-2, ISO 7176-16

GUIDE: SIZE/WEIGHT

(Recommended max user weights in kg and available sizes. Other sizes than those listed below can be ordered. Please contact your distributor)

		Width cm					
Depth cm	35	37	40	42	45	47	50
35	70		125				
40	70	100	125	125	125	150	150
42			125	125			
45	70	100	125	125	125	150	150
50		125	125	125	125	150	150