

# Curera® Consil Seat Cushion

Pressure redistribution seat cushion



Curera® Consil Seat Cushion with Standard Cover

# Curera® Consil Seat Cushion

Curera® Consil is an anatomically designed seat cushion with copyrighted design<sup>(1)</sup>. It is comfortable and pressure redistributing and provides the user with good comfort and sitting position. The anatomical design distributes the pressure away from the underside of the thighs, provides good pressure redistribution for the seat bones and reduces the risk of shearing. Side supports, leg divider and seat well provide a stable sitting position for people who need extra support.

With recesses for lumbar supports and rounded rear edge, it also provides a good fit in the wheelchair. Curera® Consil is suitable for people sitting for shorter periods. Curera® Consil may be used for wheelchairs with rigid seat base or adjustable seat base.

The cushion comes with a durable cover in grey polyester plush with friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.



Always read the instructions for use prior to use.

Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2023. All rights reserved. Last update: 2023-01-10.

## AVAILABLE COVER

### Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

### Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

### Hood cover

Made in 2-way stretch. The fabric has good flexibility and friction ability. Elastic on all sides.

### Deluxe Cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

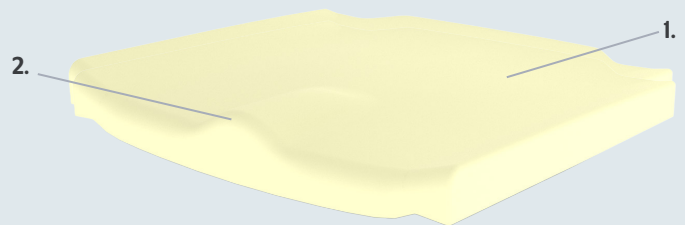

1. SEAT WELL (VISCO)

2. LEG DIVIDER

(1) Copyrighted design no 81423.

PROPERTIES	FACTS
<ul style="list-style-type: none"> <li>Anatomically moulded</li> <li>Antislip underside (Hood cover)</li> <li>Easy to clean</li> </ul>	<ul style="list-style-type: none"> <li>Height: 5 cm</li> <li>Warranty: 2 years</li> <li>Material in the cushion: Polyurethane foam</li> <li>Material in the cover: Polyester/polyurethane</li> <li>Classification: ISO-kod 181006 - Pressure redistribution</li> <li>Fire: SS-EN 1021-1, SS-EN 1021-2, ISO 7176-16</li> </ul>

## GUIDE: SIZE/WEIGHT

(Recommended max user weights in kg and available sizes. Other sizes than those listed below can be ordered. Please contact your distributor.)

Bredd cm							
Djup cm	35	37	40	42	45	47	50
35	70						
40	70	90	90	125	125	125	150
45		90	90	125	125	150	150
50					125	150	150