

Curera® Contur EX Seat Cushion

Comfortable seat cushion



Curera® Contur EX Seat Cushion with Standard Cover

Curera® Contur EX Seat Cushion

Curera® Contur EX is a comfortable and pressure redistribution seat cushion made from an anatomically moulded base of polyurethane foam with memory foam in the seat well.

The anatomical design distributes the pressure away from the underside of the thighs, provides a good and stable sitting position, good pressure redistribution for the seat bones and reduced risk of shearing. Recesses for lumbar supports, rounded rear edge and leg divider also give it good stability and fit in a wheelchair. Curera® Contur EX provides a good level of sitting comfort and is suitable for people sitting for both longer and shorter periods. Curera® Contur EX may be used for wheelchairs with rigid seat base or adjustable seat base.

Curera® Contur EX is easy to clean. The cushion comes with a durable cover in grey polyester plush with high friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.

AVAILABLE COVERS

Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

Hood cover

Made in 2-way stretch. The fabric has good flexibility and friction ability. Elastic on all sides.


Deluxe Cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

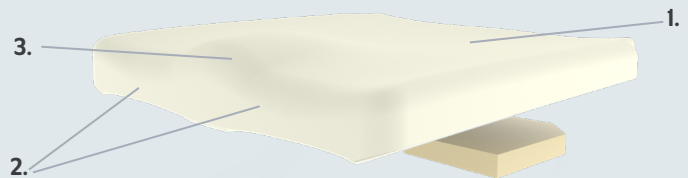

- 1. SEAT WELL
- 2. RECESSES FOR LEGS
- 3. LEG DIVIDER



PROPERTIES	FACTS
<ul style="list-style-type: none"> • Anatomically moulded • Easy to clean 	<ul style="list-style-type: none"> • Height: 5 cm • Warranty: 2 years • Material in the cushion: Polyurethane foam, memory foam • Material in the cover: Polyester/polyurethane • Classification: ISO-code 18 10 06 - Pressure redistribution • Fire: SS-EN 1021-1, SS-EN 1021-2
 Always read the instructions for use prior to use. <small>Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2023. All rights reserved. Last update: 2023-01-10.</small>	

GUIDE: SIZE/WEIGHT

(Recommended max user weights in kg and available sizes. Other sizes than those listed below can be ordered. Please contact your distributor.)

		Width cm						
Depth cm	35	37	40	42	45	47	50	
35	125		125					
40	125	125	125	125	125	125	150	
45	125	125	125	125	125	150	150	
50		125	125	125	125		150	