

Curera[®] Tremix Seat Cushion

Comfortable pressure prevention seat cushion



Curera[®] Tremix with Standard Cover

Curera® Tremix Seat Cushion

Curera® Tremix is anatomically built and used as an aid for pressure prevention. The structure of the cushion and the anatomically adapted shape with recesses for tail vertebra, splitter and lateral support make it mould to the user. It contributes to distribute the pressure on the underside of the thighs, provides a good seating position, good relief for sitting bones and tail vertebra, and helps to reduce the risk of shearing forces.

With recesses for sitting bones and rounded rear edge, it also provides a good fit in the wheelchair. The seating comfort of the Curera® Tremix is good and is suitable for people who sit both longer and shorter time. For extra good seating function it is equipped with reinforced bottom plate.

Curera® Tremix has a liquid-tight finish that does not allow liquids and is easy to clean. The cushion follows a durable coating of grey polyesterplysch with high friction. All Curera® seat cushions and covers do not contain brominated flame retardants or latex.



Always read the instructions for use prior to use.

Care of Sweden is a registered trademark of care of Sweden AB. © Care of Sweden AB, 2021. reiterated. Last Updated: 2021-06-07.

AVAILABLE COVERS

Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

Hood cover

Made in 2-way stretch. The fabric has good flexibility and friction ability. Elastic on all sides.

Deluxe Cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

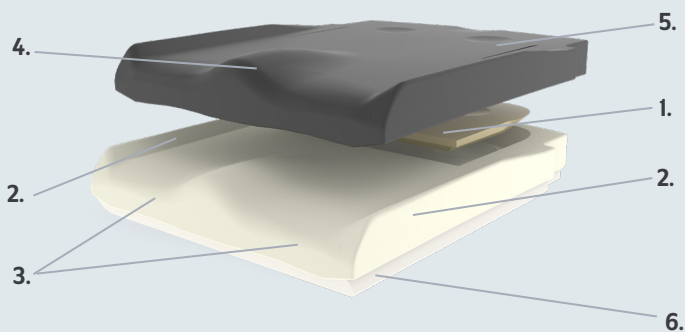

- | | |
|----------------------|----------------------------|
| 1. SEAT WELL | 4. LEG DIVIDER |
| 2. LATERAL SUPPORT | 5. WATERPROOF SURFACE |
| 3. RECESSES FOR LEGS | 6. REINFORCED BOTTOM PLATE |

PROPERTIES	FACTS
<ul style="list-style-type: none"> Anatomically built Moulds after the user Pressure preventing Easy to clean 	<ul style="list-style-type: none"> Height: 7 cm Warrant: 2 years Material in Cushion: Polyurethane foam, Visco foam, vulcanized rubber Material in cover: Polyester/polyuretan Classification: ISO Code 04 33 03-Pressure relief Fire: Complies with ISO 7176-16

GUIDE: SIZE/WEIGHT

(Recommended maximum user weights in kg and available sizes. Sizes other than those listed below are available to order. Please contact your distributor.)

Width cm		35	37	40	42	45	47	50	52	55
Depth cm	35	70								
40	70	105	115	130	130	130	150	150	160	
45		105	115	130	130	150	150	150	160	
50			115	130	130	150	150	150	160	