

CuroCell® Nova

Simple & functional mattress overlay system




CuroCell Nova with cover Olivia

FACTS	OPTIONAL
<ul style="list-style-type: none"> • Type: Dynamic mattress overlay system • Height: 13 cm • Function: Prevention and treatment up to and including category 2⁽¹⁾ • Recommended user weight: 0-150 kg • Input voltage: 230 V/50 Hz • Sound level, pump, max: 21 dBA • Warranty: 2 years • Machine wash max 95 °C, tumble drying • Cleaning of cover: Wipe with cleaning agent and/or disinfectants 	<ul style="list-style-type: none"> • Hygienic cover with welded seams • Durable bag for transportation and storage
	<p>SIZES</p> <p>80x200 cm.</p>

CuroCell® Nova

The CuroCell® Nova is used as an aid in the prevention and treatment of pressure ulcers up to category 2⁽¹⁾. It is a mattress overlay system that is intended for community and home care and is easy to handle and operate. The one-piece mattress has air cells of polyurethane and a built-in heel function.

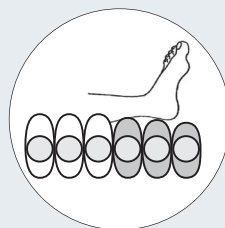


REFERENCES

(1) European Pressure Ulcer Advisory Panel, National Pressure Injury Advisory Panel and Pan Pacific Pressure Injury Alliance. Prevention and Treatment of Pressure Ulcers/Injuries: Clinical Practice Guideline. The International Guideline. Emily Haesler (Ed). EPUAP/NPIAP/PPPIA: 2019.

(2) CDC (2003). Guidelines for Environmental Infection Control in Health-Care facilities, updated 2019.



FEATURES CONTROL UNIT	FEATURES MATTRESS
<ul style="list-style-type: none"> • Silent running and durable pump • Energy efficient • Static function and Max firm function • Comfort settings on dial, 0-150 kg 	<ul style="list-style-type: none"> • CPR function – rapid deflation • Soft, elastic, vapour permeable⁽⁴⁾ hygienic cover • Good cleaning properties • Easy to handle and maintain • Elastic fastening band <p> Always read the instruction for use prior to usage of the product.</p> <p>Care of Sweden and CuroCell are trademarks of Care of Sweden AB. © Care of Sweden AB, 2021. All rights reserved. Last update: 2021-06-08.</p>