


Optimal 5zon[®] Plus

Preventative mattress replacement for plus-size use

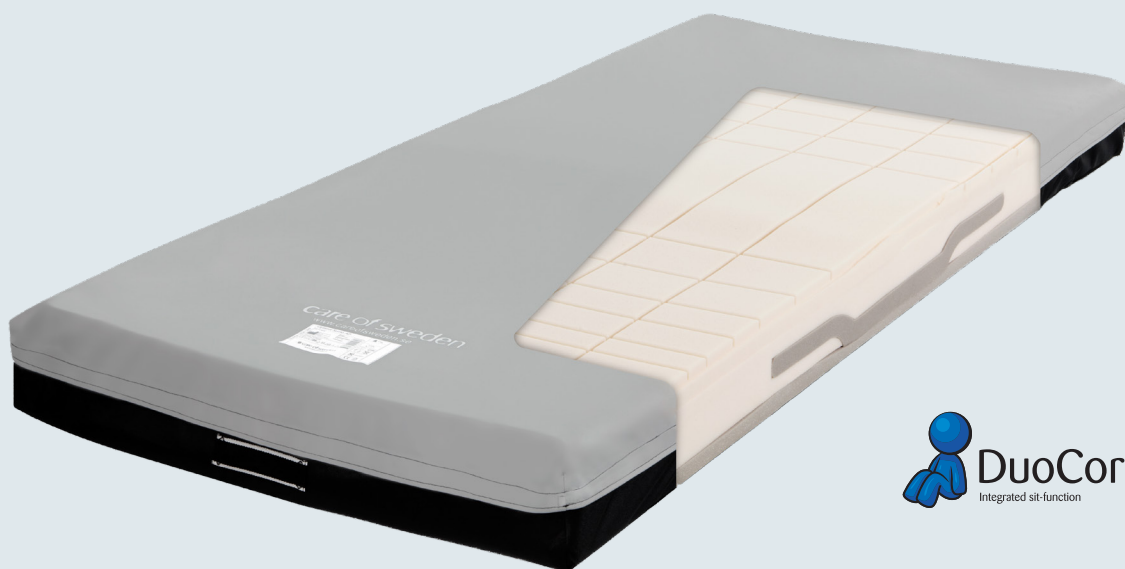

Optimal 5zon Plus with top part Olivia/bottom part 5zon Plus

Optimal 5zon[®] Plus

Optimal 5zon[®] Plus is used as an aid in the prevention and treatment of pressure ulcers up to category 2⁽¹⁾, and is intended for all healthcare environments. It is a mattress replacement that has been specially designed for plus-size use.

It is equipped with our DuoCore[®] sit function, which provides assistance and relief when seated, as well as good comfort in a lying position. In addition, the DuoCore[®] is also beneficial when using a sliding board or other relocation aids.

The cover is equipped with a divisible, drip-shield zip and functional handles for simple transportation of the mattress.

Care of Sweden's mattresses are supplied with a four way stretch hygienic outer cover which conforms to relevant standards for hygiene, quality and safety⁽²⁾⁽⁵⁾.



All covers are manufactured from modern, recyclable materials and are vapour permeable⁽³⁾, which contributes to reduce the risk of skin maceration.

SIZE

85/90/105/120x200/210 cm.

If you need any other sizes please contact your distributor.

REFERENCES

(1) European Pressure Ulcer Advisory Panel, National Pressure Injury Advisory Panel and Pan Pacific Pressure Injury Alliance. Prevention and Treatment of Pressure Ulcers/Injuries: Clinical Practice Guideline. The International Guideline. Emily Haesler (Ed). EPUAP/NPIAP/PPPIA: 2019

(2) CDC (2003), Guidelines for Environmental Infection Control in Health-Care facilities, updated 2019.

(3) SS-EN ISO 15496:2004, DIN 53122-1.

(4) Patent no. 008522741.

(5) SS 876 00 20:2017, Healthcare textiles - Mattresses - Specifications and requirements



Always read the instruction for use prior to usage of the product.

Care of Sweden, Optimal 5zon and DuoCore are trademarks of Care of Sweden AB. © Care of Sweden AB, 2021. All rights reserved. Last update: 2021-06-02.

FEATURES & CHARACTERISTICS	FACTS
<p>User weight For plus-size use, recommended user weight is up to 300 kg.</p> <p>Pressure Elastic material and a zone and module system gives good pressure equalisation.</p> <p>Shear forces The mattress is specially designed to reduce the effects of shear forces.</p> <p>Sit function Our patented⁽⁴⁾ sit function DuoCore[®] is designed to give support when sitting and also to provide good comfort when lying.</p> <p>Moisture The hygienic cover is vapour permeable⁽³⁾, which contributes to reduce the risk of skin maceration.</p> <p>Hygiene The cover is easy to clean and has a good hygienic standard⁽²⁾⁽⁵⁾.</p>	<ul style="list-style-type: none">• Core: HR-foam 65+50 kg/m³• Height: 15 cm• Function: an aid for prevention and treatment up to and including category 2⁽¹⁾• Recommended user weight: up to 300 kg (Optimal 5zon[®] is recommended for user weights up to 180kg).• Fire: EN 597-1, -2• Warranty: 2 years• Cleaning of cover: wipe with cleaning agent and/or disinfectants.• Machine wash max 95 °C, tumble drying