

Curera[®] Consept Latus Seat Cushion

Seat cushion for heavy users



Curera[®] Consept Latus with Standard Cover

Curera® Consept Latus Seat Cushion

Curera® Consept Latus is a comfortable and pressure redistribution seat cushion intended for heavy users, made from an anatomically moulded base of polyurethane foam with memory foam filling on the top side. Fitted with reinforced bottom plate, side supports, leg divider and seat well, the cushion provides a stable seat for people who need extra support.

The anatomical and adapted design distributes the pressure away from the underside of the thighs, provides a good sitting position, good pressure redistribution for the seat bones and caudal vertebra and reduced risk of shearing. Recesses for lumbar supports and rounded rear edge also give it good stability and fit in a wheelchair. It is a comfortable seat cushion and is suitable for people sitting for both longer and shorter periods.

Concept Latus has a waterproof surface that prevents liquids from passing through and is easy to clean. The cushion comes with a durable cover in grey polyester plush with high friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.



Always read the instructions for use prior to use.

Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2017. All rights reserved. Last update: 2021-05-31.

AVAILABLE COVERS

Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

Deluxe Cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

Onyx Cover

Material made of 4-way stretch on the top side and an anti-slip underside. Comes with handles on the underside. Available with a mesh layer for a better microclimate and a generally higher comfort.

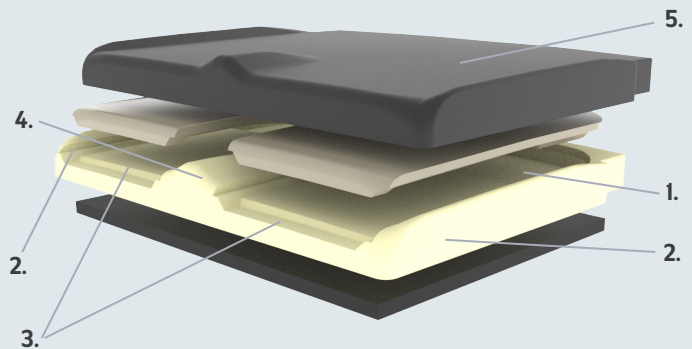

- 1. SEAT WELL
- 2. LATERAL SUPPORTS (SIDE SUPPORTS)
- 3. RECESSES FOR LEGS
- 4. LEG DIVIDER
- 5. WATERPROOF SURFACE

PROPERTIES	FACTS
<ul style="list-style-type: none"> • Anatomically moulded • Easy to clean • Bottom plate of vulcanised rubber 	<ul style="list-style-type: none"> • Height: 7 cm • Warranty: 2 years • Material in the cushion: Polyurethane foam, memory foam, vulcanised rubber • Material in the cover: Polyester/polyurethane • Classification: ISO-kod 181006 - Pressure redistribution • Fire: SS-EN 1021-1, SS-EN 1021-2, ISO 7176-16

GUIDE: SIZE/WEIGHT

(Recommended max user weights in kg and available sizes. Other sizes than those listed below can be ordered. Please contact your distributor.)

Bredd cm					
Djup cm	50	52	55	57	60
40		250	250	250	370
42	250				
45	250	250	250		370
47					370
50	250	250	250	370	370
55	250		250	370	370