

Curera® Contur Seat Cushion

Comfortable and popular seat cushion

Curera® Contur is a comfortable and pressure redistribution seat cushion made from an anatomically moulded base of polyurethane foam with memory foam in the seat well.

The anatomical design distributes the pressure away from the underside of the thighs, provides a good and stable sitting position, good pressure redistribution for the seat bones and reduced risk of shearing. Recesses for lumbar supports, rounded rear edge and leg divider also give it good stability and fit in a wheelchair. Contur provides a good level of sitting comfort and is suitable for people sitting for both longer and shorter periods. Curera® Contur may be used for wheelchairs with rigid seat base or adjustable seat base.

Contur has a waterproof surface that prevents liquids from passing through and is easy to clean. The cushion comes with a durable cover in grey polyester plush with high friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.



Curera® Contur Seat Cushion with standard cover

AVAILABLE COVERS

Standard cover


Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic cover

A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

Deluxe cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

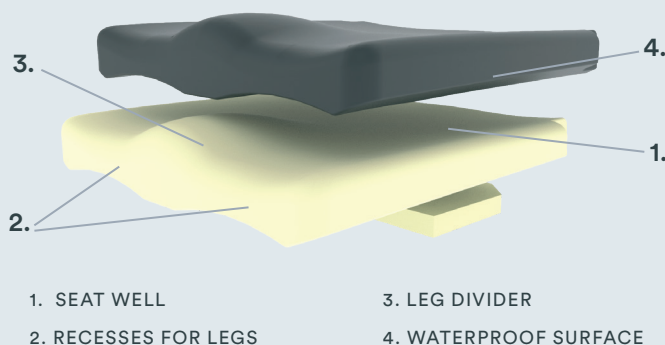

Always read the instructions for use prior to use.

			PROPERTIES	FACTS
  			<ul style="list-style-type: none"> Anatomically moulded Easy to clean Accessories: Levelling plate 	<ul style="list-style-type: none"> Height: 5 cm Warranty: 2 years Material in the cushion: Polyurethane foam, memory foam, vulcanised rubber Material in the cover: Polyester/polyurethane Classification: ISO-code 122506-03 - Pressure redistribution Fire: Complies with EN 1021-1, EN 1021-2, ISO 7176-16
<small>Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2025. All rights reserved. Last update: 2025-04-14.</small>				

PASSFORMSGUIDE:

(Recommended max user weights in KG and available standard sizes. Other sizes than listed below can be ordered. Please contact your distributor.)

Width cm											
Depth cm	30	33	35	38	40	42	45	47	50	52	
30	125	125									
35	125	125	125		125						
40	125	125	125	125	125	125	125	125	150	150	
42					125	125	125	125	150		
45			125	125	125	125	125	150	150	150	
47					125	125	125				
50				125	125	125	125		150	150	
55									150		

Please note!
Make sure to measure where the arrows indicate.

