

# Curera® Hemi Seat Cushion

## Seat cushion with special design

Curera® Hemi is a pressurised, specially designed cushion with registered design<sup>(1)</sup> and suitable for people who need to be able to kick their way using the wheelchair. It is constructed of an anatomically molded base with inserts of Visco foam from the top. For better mobility, the cushion is at the front and front half submerged on one or both sides. Equipped with lateral support, leg dividers and seat well, the cushion provides a stable seat for people who need extra support.

The anatomical and custom shape distributes the pressure on the underside of the thighs, provides a good seating position, good relief for sitting bones and tail vertebra, and reduced risk of shearing forces. With recesses for lumbar supports and rounded rear edge, it also provides a good fit in the wheelchair. The seating comfort of the Hemi is high and is suitable for people who sit both longer and shorter time.

It has a liquid-proof finish that does not let through liquids and is easy to clean. The cushion follows a durable coating of grey polyesterplysch with high friction.

Curera® Hemi is available in three variants, submerged on the right or left side or on both sides. All Curera® seat cushions and covers do not contain brominated flame retardants or latex.



Curera® Hemi Seat Cushion with Standard Cover

## AVAILABLE COVERS

### Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

### Hygienic Cover


A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

### Hood cover

Made in 2-way stretch. The fabric has good flexibility and friction ability. Elastic on all sides.

### Deluxe Cover

Plush polyester on the outside and a self-ventilating material inside. Elastic and flexible material on the top side with an anti-slip underside.

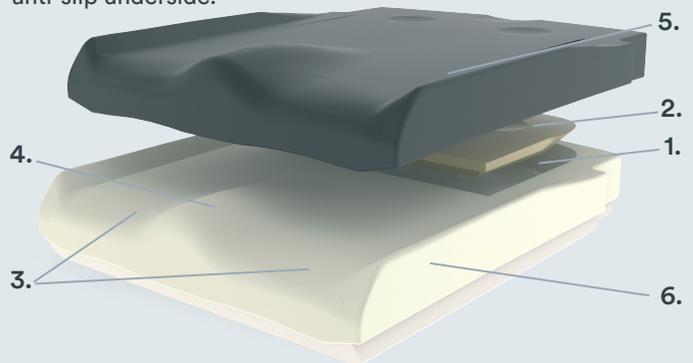

- |                               |                        |
|-------------------------------|------------------------|
| 1. SEAT WELL                  | 4. BONE SPLITTER       |
| 2. RECESSES FOR TAIL VERTEBRA | 5. FLUID-PROOF COATING |
| 3. RECESSES FOR LEGS          | 6. SIDE SUPPORT        |



Always read the instructions for use prior to use.

    <small>(1) Design protection no 79436. Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2026. All rights reserved. Last update: 2026-01-23.</small>	PROPERTIES	FACTS
	<ul style="list-style-type: none"> <li>Anatomically moulded</li> <li>Front half of the cushion lowered on one or two sides</li> <li>Easy to clean</li> </ul>	<ul style="list-style-type: none"> <li>Height: 5 cm</li> <li>Warranty: 2 years</li> <li>Material in Cushion: Polyurethane foam, Visco foam, vulcanized rubber</li> <li>Material in cover: Polyester/polyuretan</li> <li>Classification: ISO Code 122506-03 - Pressure redistribution</li> <li>Fire: SS-EN 1021-1, SS-EN 1021-2, ISO 7176-16</li> </ul>

## GUIDE: SIZE/WEIGHT

(Recommended maximum user weights in kg and available sizes. Sizes other than those described below can be ordered. Please contact your distributor.)

Depth cm	Width cm						
	35	37	40	42	45	47	50
35	70		70	115			
40	70	70	115	130	130	130	150
45		70	115	130	130	130	150
50				130	130		150

### Please note!

Make sure to measure where the arrows indicate.

