

Curera® Plantur Seat Cushion

Seat Cushion with flat top side

Curera® Plantur is a pressure redistribution seat, made from moulded polyurethane foam with a flat top side and with memory foam filling on the underside.

The seat well and recesses for legs distribute the pressure away from the underside of the thighs, providing a good sitting position, good pressure redistribution for the seat bones and reduced risk of shearing. Recesses for lumbar supports and rounded rear edge also give it good stability and fit in a wheelchair. Plantur provides comfort suitable for people sitting for both longer and shorter periods.

Curera® Plantur has a waterproof surface that prevents liquids from passing through and is easy to clean. The cushion comes with a durable cover in grey polyester plush with high friction. None of Curera® seat cushions and covers contains any bromide-based flame retardants or latex.



Curera® Plantur Seat Cushion with Standard Cover


AVAILABLE COVERS

Standard Cover

Flexible and elastic material on the top side that prevents the user sliding with an anti-slip underside.

Hygienic Cover

A liquid proof and vapour permeable cover that is elastic and flexible with an anti-slip underside.

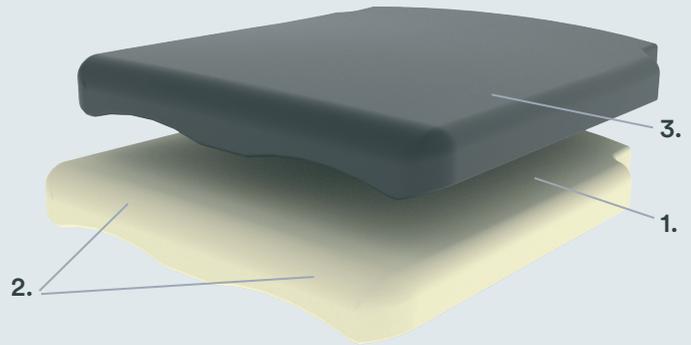

1. SEAT WELL
2. RECESSES FOR LEGS
3. WATERPROOF SURFACE



Always read the instructions for use prior to use.

|     <small>Care of Sweden is a trademark of Care of Sweden AB. © Care of Sweden AB, 2026. All rights reserved. Last update: 2026-01-23.</small> | PROPERTIES | FACTS |
|---|---|--|
| | <ul style="list-style-type: none"> • Easy to clean • Flat surface | <ul style="list-style-type: none"> • Height: 5 cm • Warranty: 2 years • Material in the cushion: Polyurethane foam, memory foam, vulcanised rubber • Material in the cover: Polyester/polyurethane • Classification: ISO-kod 122506-03 - Pressure redistribution • Fire: SS-EN 1021-1, SS-EN 1021-2, ISO 7176-16 |

GUIDE: SIZE/WEIGHT

(Recommended max user weights in kg and available sizes. Other sizes than those listed below can be ordered. Please contact your distributor)

| | Width cm | | | | | | |
|----------|----------|-----|-----|-----|-----|-----|-----|
| Depth cm | 35 | 37 | 40 | 42 | 45 | 47 | 50 |
| 35 | 70 | | 125 | | | | |
| 40 | 70 | 100 | 125 | 125 | 125 | 150 | 150 |
| 42 | | | 125 | 125 | | | |
| 45 | 70 | 100 | 125 | 125 | 125 | 150 | 150 |
| 50 | | 125 | 125 | 125 | 125 | 150 | 150 |

Please note!

Make sure to measure where the arrows indicate.

